



somté[®]PSG

Full PSG... absolutely anywhere

 COMPUMEDICS[®]

Whether you are looking for a recorder for attended sleep studies, or you wish to perform full polysomnography in an unattended, ambulatory setting, or both,



Full PSG for in lab convenience


- Type I Attended in-lab recordings of the full complement of PSG signals using recommended technology.
- Maximum patient freedom and comfort - The patient is not tethered to any bedside unit or computer, instead data is transferred via Bluetooth wireless transmission from the patient worn unit to the bedside station.
- Do everything you can usually do with an in-lab study – monitor patient data, record synchronized digital video, add annotation, run automatic or manual on-line data analysis. Technical reports can be ready just minutes after the recording ends.
- The bedside station provides Ethernet connectivity for remote LAN/WAN/Internet viewing and acquisition.
- The bedside station allows up to 8 external DC inputs for other channel capabilities such as pH or capnography.



Pictured actual size

somté[®] PSG

is the solution
you're looking for - the simplest and most convenient way to
meet requirements for **full PSG ... absolutely anywhere.**

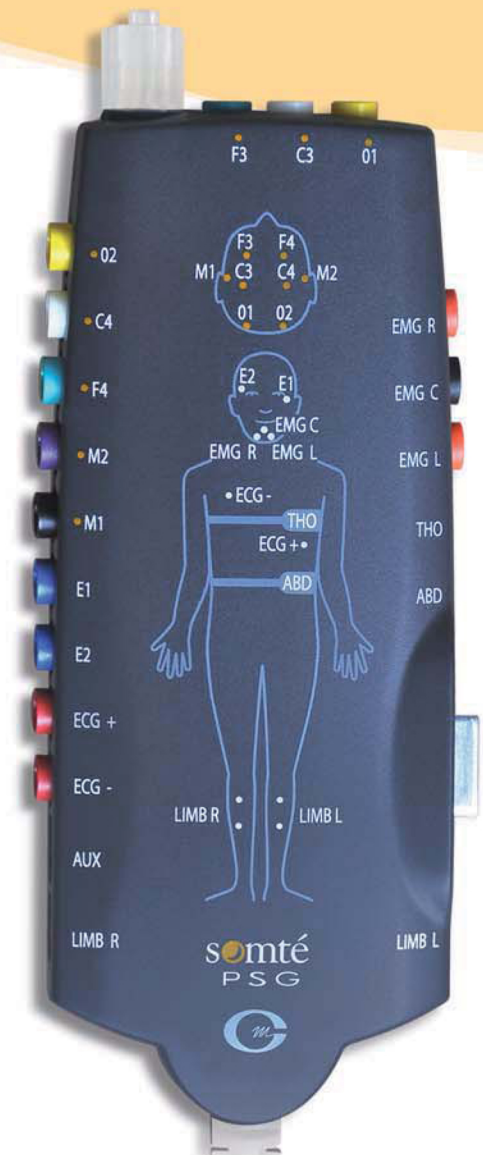

Full PSG for ambulatory recording

- The simple patient-worn unit allows rapid setup and use of the system anywhere.
- Fully self contained. No PC or laptop is required to setup and acquire data, although data cards can be pre-programmed to automatically start and stop recording at pre-set times, and add patient identification to the study.
- Integrated LCD screen to check the quality of signals before or during recording.
- Unique signal status display for instant signal integrity and impedance checks. Signal integrity can be displayed during data review for quality assurance.
- Integrated storage card (Compact Flash) for up to 24 hours data collection.
- Uses analogue differential amplifiers for the simplest and most robust recording of physiological signals, even in less than ideal situations.

Technical Specifications

25 Available channels

- 6 x EEG
- 2 x EOG
- 2 x Chin EMG
- ECG – with Thoracic Impedance option
- Airflow from Pressure transducer
- Snore
- Pressure (for CPAP/ Bilevel)
- Thoracic Inductive RIP band
- Abdominal Inductive RIP band
- 2 x Limb EMG
- Position
- Oximeter provides SpO₂, Pulse Rate, Plethysmography waveform and Oximeter probe status
- Aux selectable as thermistor, external position sensor, Aux AC or Aux DC

Weight

- 205gm (7.2 ounces) (excluding batteries)
- **Recorder Size:** Width 113mm x Height 65mm x Depth 30mm (including belt clip), (4.4 x 2.5 x 1.2 inches)
- **Patient Input Box:** Width 53mm x Height 133mm x Depth 25mm (including belt clip), (2.1 x 5.2 x 1 inches)

Flexible Signal Options

- Use a cannula and pressure transducer and/or thermistor to record airflow
- Use Inductive Plethysmography or thoracic impedance technology to record respiratory effort
- Use position sensor integrated into headbox or an external position sensor if preferred

Sampling Rate

- 1-200Hz

Wireless Technology

- Uses Bluetooth Wireless technology (Class 2) for patient freedom
- Up to 10m range from bedside station to patient - dependant on environment

Power

- 2 x AA alkaline or NiMH batteries

Recording Time

- For up to 24 hours of recording

Event Button

- Any button press during recording is recorded as an event

Electrode/Sensor status

- Electrode/sensor status indicator before recording

Compatible Software

- **Data Card Manager** allows selection of options, pre-set start and stop recording times, set up patient information
- **PSG Online** for monitoring attended studies via Bluetooth. Not required for unattended studies
- **ProFusion PSG** for data analysis
- **Digital Video** can be used with PSG Online
- **ProFusion neXus Lab Manager** for data and workflow management

Somté PSG System

Part Number: 9023-0002-01



Optional EZ Vest
for use with Somté PSG

Since 1987, Compumedics Limited has emerged as a world leader in the design, manufacture, and sales of diagnostic sleep, neurophysiology, research and transcranial Doppler systems. Compumedics design expertise in the area of portable physiologic recorders and wireless data transmission is just one reason for our continued growth. Compumedics is improving people's lives through innovative products, exceptional service and a full range of consumables and accessories brought to you by dedicated people who care.

Compumedics divisions:



For more information please contact:

Compumedics Limited, Australia:
30-40 Flockhart Street
Abbotsford VIC 3067
Australia
Tel: +61 3 8420 7300
Fax: +61 3 8420 7399
Free Call: 1800 651 751

Compumedics USA, Inc.:
6605 West WT Harris Blvd, Suite F
Charlotte, NC 28269
Toll Free: +1 877 717 3975
Ph: +1 704 749 3200
Fax: +1 704 749 3299

Compumedics Germany GmbH:
Global HQ for Compumedics DWL
Josef-Schüttler-Strasse 2
D-78224 Singen, Germany
Ph: +49 7731 79 76 9-0
Fax: +49 7731 79 76 9-99

Your local Compumedics Distributor:

All specifications are subject to change without notice. Please contact your Compumedics representative for latest technical information, pricing and product availability. Compumedics Somté PSG, COMPUMEDICS and the Compumedics logo are all registered trademarks of Compumedics.

 **COMPUMEDICS®**
www.compumedics.com